



Compressing Machine for Chamber Filters “ FilTrex ”

FEATUERS:

- The compressing machine was designed to operate with our FilTrex Chamber Filters.
- It can compatible with other Syringe-less Filters or Vial Filters too.
- 10 vials can do simultaneously at once compressing , resulting in a fast sample preparation.
- The material of body made of Stainless steel by laser cutting and the toggle head made of cast iron.
- The rack of vial can move in/out from the machine. so it can move of all 10 samples to analytical procedure ,instrument or autosampler.
- Each hole of vial can see the liquid inside the vial pass through the window view.
- It is easy to operate by compress the hand of toggle head down for the filtration process.

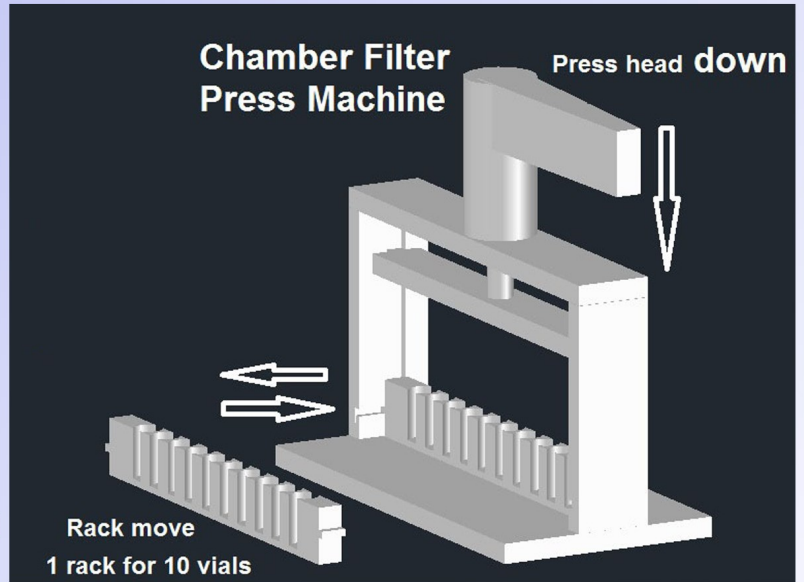

Description	Part no.	Price
Compressing Machine for Chamber Filters ,10 positions, stainless steel body , 4 cm. stoke ,weight 2.0 Kg.	FTC-CM010-001	

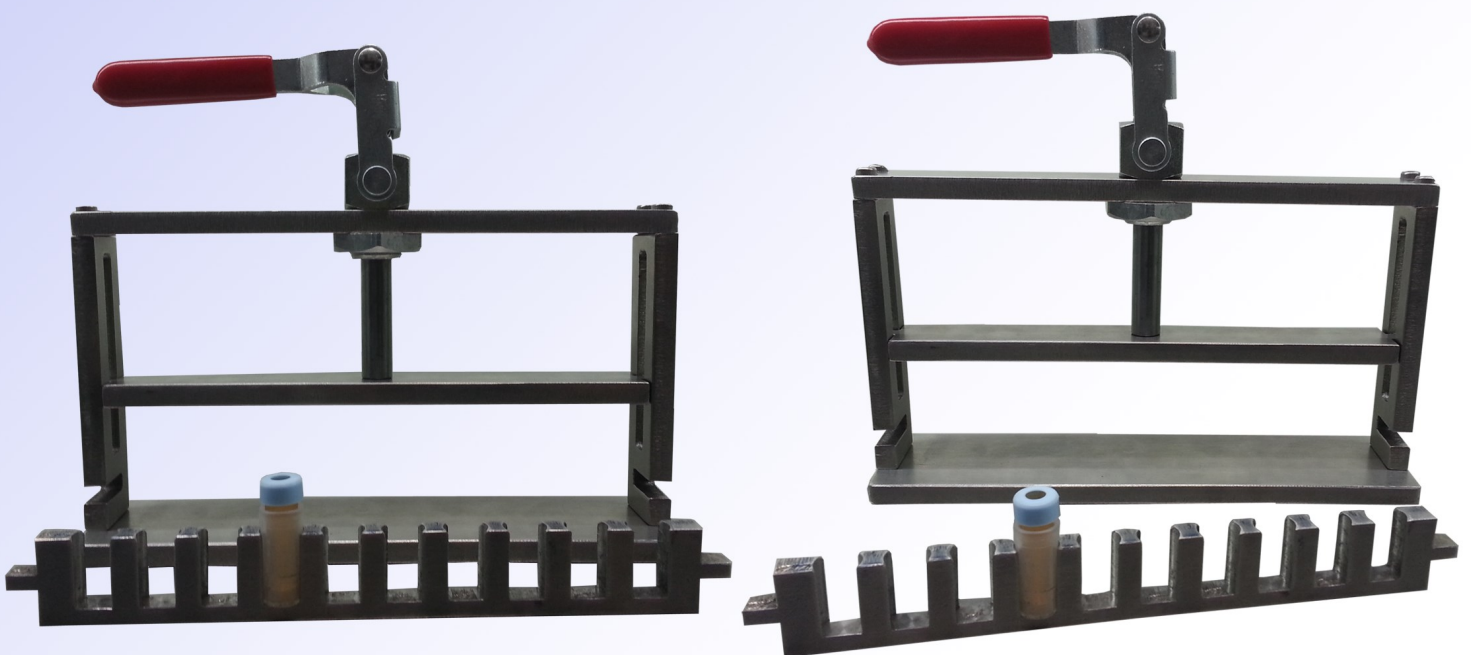

NOTE : We reserve the right to alter the specifications of the products described in this leaflet without prior notice.
www.lab-preparation.com